

Release Note

Product Name: PM400

Version: V2.2

Release Date: February, 2022

What's New?

This firmware version release implements enhancements in the PM400 monitors with the following changes:

- Minor improvements for manufacturing
- Minor bug fixes


Installation Instructions

Installation

This firmware version requires PM Link software version 1.7.0 or above.

Installing the new version:

1. Download the new firmware .bin file (Make sure to remember where the new firmware is saved on the PC)
2. Power up the PM400
3. Connect the PM Link to the PC
4. Place the PM400 in front of the PM Link aligning the IR window on the PM400 with the PM Link
5. Launch the PM Link software
6. Click the OPEN button to establish communication with the PM400 and the PM Link PC Read Complete popup window will appear. Click the OK button.

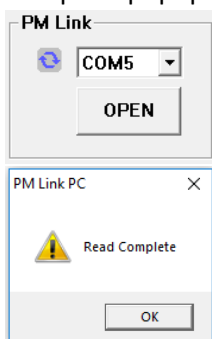


In the Basic Operation section, the current firmware version will be displayed.

Firmware Version	1.7
------------------	-----

7. Click on Firm icon and the Firmware Upgrade popup window will appear




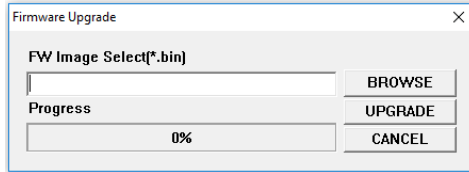


8. Click on the BROWSE button and locate the new firmware .bin file (Typically in Downloads)
9. Click on the appropriate .bin file to select it, then click the Open button
10. Ensure that the PM400 is still properly located in front of the PM Link
11. Click on the UPGRADE button




12. The display on the PM400 will display IAP and the Progress bar on the software will indicate the progress of the firmware installation. (It should take less than 30 seconds to update the firmware)

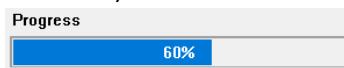



Once the firmware has upgraded, the PM400 will automatically reboot and the Firmware Upgrade Success popup window will be displayed on the PC. Click the OK button.

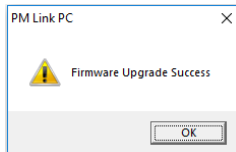


13. Click the X button on the Firmware Upgrade popup window to close the upgrade window

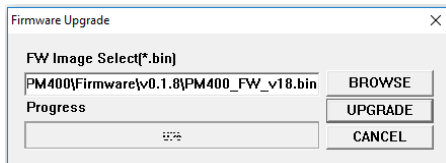

14. Click the OPEN button (as described in step 6) to establish communication with the PM400 and the PM Link PC Read Complete popup window will appear. Click the OK button.

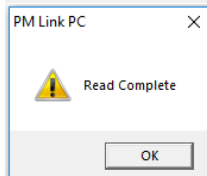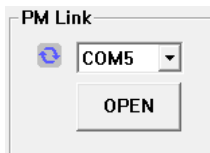

In the Basic Operation section, verify that the new firmware has been loaded.

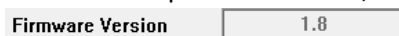

15. The firmware upgrade is complete

Contact Technical Support with any questions

Phone: (844) 325-3050

Email: support@aimsafety.com